

Exhibit 4

- (54) Title
Recombinant super-compound interferon
- (51)⁶ International Patent Classification(s)
A61K 038/21 A61P 031/12
C12N 015/63 A61P 001/16
C12N 015/20 C12N 015/70
- (21) Application No: **2003248419** (22) Application Date: **2003.09.26**
- (43) Publication Date : **2003.11.06**
- (43) Publication Journal Date : **2003.11.06**
- (62) Divisional of:
2002235707
- (71) Applicant(s)
Sichuan Biotechnology Research Center
- (72) Inventor(s)
Zhang, Renhui; Guo, Rongbing; Wei, Guangwen⁷
- (74) Agent/Attorney
Patent Attorney Services, 26 Ellingworth Parade, BOX HILL, VIC, 3128

AUSTRALIA

Patents Act 1990

**COMPLETE SPECIFICATION
FOR A STANDARD PATENT
ORIGINAL**

Applicant(s): SICHUAN BIOTECHNOLOGY RESEARCH CENTER

Actual Inventor(s): Guangwen WEI, Rongbing GUO and Renhual ZHANG

Address for Service: **PATENT ATTORNEY SERVICES**
26 Ellingworth Parade
Box Hill Victoria 3128
Australia

Title: RECOMBINANT SUPER-COMPOUND INTERFERON

Associated Provisional Applications: No(s):

The following statement is a full description of this invention, including the best method of performing it known to me/us:-

1

RECOMBINANT SUPER-COMPOUND INTERFERON

The application is a continuation-in-part application of International Patent Application No. PCT/CN02/00129, filed on 28 February 2002, which claims priority of Chinese Application No. 01104367.9, filed on 28 February 2001, the contents of which are incorporated by reference here into this application.

Throughout this application, various references are referred to. Disclosures of these publications in their entireties are hereby incorporated by reference into this application in order to more fully describe the state of the art to which this invention pertains.

FIELD OF THE INVENTION

This invention is related to a recombinant super-compound interferon (rSIFN-co) with changed spatial configuration. One characteristic of rSIFN-co in this invention is that it cannot only inhibit DNA (deoxyribonucleic acid) duplication of the Hepatitis B virus but also the secretion of HBsAg and HBeAg.

BACKGROUND OF THE INVENTION

rSIFN-co is a new interferon molecule constructed with the most popular conservative amino acid found in natural human α -IFN subtypes using genetic engineering methods. United States Patent Nos. 4,695,623 and 4,897,471 have described it. rSIFN-co had been proved to have broad-spectrum IFN activity and virus- and tumor-inhibition and natural killer cell activity. United States Patent No. 5,372,808 by Amgen, Inc. addresses treatment rSIFN-co. Chinese Patent No. 97193506.8 by Amgen, Inc. addresses re-treatment of rSIFN-co on Hepatitis C. Chinese Patent No. 98114663.5 by Shenzhen Jiusheng Bio-engineering Ltd. addresses treatment

of rSIFN-co on Hepatitis B and Hepatitis C.

The United States Food and Drug Administration (FDA) authorized Amgen to produce rSIFN-co with E. Coli. for
5 clinical Hepatitis C treatment at the end of 1997.

Hepatitis B patients can be identified when detecting HBsAg and the HBeAg. α -IFN is commonly used in clinics to treat Hepatitis B. IFN binds superficial cell membrane receptors,
10 inhibiting DNA and RNA (ribonucleic acid) duplication, including inducing some enzymes to prevent duplication of the virus in hepatitis-infected cells. All IFNs can inhibit only the DNA duplication of viruses, not the e and s antigen.

15

This disclosure describes recombinant super-compound interferon, method to produce the same and uses thereof.

The above references to and descriptions of prior proposals
20 or products are not intended to be, and are not to be construed as, statements or admissions of common general knowledge in the art in Australia.

SUMMARY OF THE INVENTION

This invention provides a recombinant super-compound interferon or an equivalent thereof with changed spatial configuration. An equivalent is a molecule which is similar in function to the super-compound interferon. The super-compound interferon possesses anti-viral or anti-tumor activity. This invention also provides an artificial gene codes for the super-compound interferon or its equivalent.

This invention provides a process for production of recombinant super-compound interferon comprising introducing an artificial gene with selected codon preference into an appropriate host, culturing said introduced host in an appropriate condition permitting expression of said super-compound interferon and harvesting the expressed super-compound interferon.

This invention provides a composition comprising the recombinant super-compound interferon or its equivalent and a suitable carrier. This invention further provides a pharmaceutical composition comprising the recombinant super-compound interferon or its equivalent and a pharmaceutically acceptable carrier.

This invention provides a method for treating viral diseases or tumor in a subject comprising administering to the subject an effective amount of the super-compound interferon or its equivalent.

This invention provides the above-described method wherein super-compound interferon was administered via oral, vein injection, muscle injection, peritoneal injection, subcutaneous injection, nasal, mucosal administration, by inhalation via an inspirator.

DETAILED DESCRIPTION OF THE FIGURES

5 **Figure 1.** rSIFN-co cDNA sequence designed according to E. Coli. codon usage and deduced rSIFN-co amino acid sequence

Figure 2. Sequence of another super-compound interferon

10 **Figure 3.** Diagram of pLac T7 cloning vector plasmid

Figure 4. Diagram of pHY-4 expression vector plasmid

Figure 5. Construction process of expression plasmid pHY-5

15 **Figure 6-A.** Circular Dichroism spectrum of Infergen®

Spectrum range: 250nm - 190nm

Sensitivity: 2 m°/cm

Light path: 0.20 cm

20 Equipment: Circular Dichroism J-500C

Samples: contains 30µg/ml IFN-con1, 5.9 mg/ml of NaCl and 3.8 mg/ml of Na₂PO₄, pH7.0.

INFERGEN® (interferon alfacon-1), made by Amgen Inc., also
25 known as consensus interferon, is marketed for the treatment of adults with chronic Hepatitis C virus (HCV) infections. It is currently the only FDA approved, bio-optimized interferon developed through rational drug design and the only interferon with data in the label specifically
30 for non-responding or refractory patients. InterMune's sales force re-launched Infergen® in January 2002 with an active campaign to educate U.S. hepatologists about the safe and appropriate use of Infergen®, which represents new hope for the more than 50 percent of HCV patients who fail
35 other currently available therapies. See

<http://www.intermune.com/vt/itmn/infergen>, 8/27/2003

Figure 6-B. Circular Dichroism spectrum of Infergen® From
Reference [Journal of Interferon and Cytokine Research.
16:489-499(1996)]

Figure 6-C. Circular Dichroism spectrum of rSIFN-co

Spectrum range: 320nm-250nm

Sensitivity: 2 m°/cm

Light path: 2cm

Equipment: Circular Dichroism J-500C

Sample: contains 0.5mg/ml rSIFN-co, 5.9 mg/ml of NaCl and
3.8 mg/ml of Na₂PO₄, pH7.0.

Figure 6-D. Circular Dichroism spectrum of rSIFN-co

Spectrum range: 250nm - 190nm

Sensitivity: 2 m°/cm

Light path: 0.20 cm

Equipment: Circular Dichroism J-500C

Sample: contains 30pg/ml rSIFN-co, 5.9 mg/ml of NaCl and
3.8 mg/ml of Na₂PO₄, pH7.0.

Clearly, as evidenced by the above spectra, the secondary
or even tertiary structure of rSIFN-co is different from
Infergen®.

DETAILED DESCRIPTION OF THE INVENTION

This invention provides a recombinant super-compound interferon or an equivalent thereof with changed spatial configuration. This invention reveals that protein with same primary sequence might have different biological activities. As illustrated in the following example, this invention disclosed two proteins with identical amino acid sequence but with different activities. This activity may sometimes become improved efficacy and sometimes, the protein with changed spatial configuration would reveal new function.

An equivalent is a molecule which is similar in function to the compound interferon. An equivalent could be a deletion, substitution, or replacement mutant of the original sequence. Alternatively, it is also the intention of this invention to cover mimics of the recombinant super-compound interferon. Mimics could be a peptide, polypeptide or a small chemical entity.

The interferon described herein includes but is not limited to interferon α , β , or ω . In an embodiment, it is IFN-1 α , IFN-2b or other mutants.

In an embodiment, the super-compound interferon disclosed has higher efficacy than the interferon described in U.S. Patent Nos. 4,695,623 or 4,897,471. This super-compound interferon is believed to have unique secondary or tertiary structure. (See e.g. Figure 6)

The super-compound interferon described herein has special structure change(s) resulting from the changes of its production process.

The above-described super-compound interferon may be

produced by a high efficiency expression system which uses a special promoter. In an embodiment, the promoter is P_{lac}. As it could be easily appreciated by other ordinary skilled artisan, other inducible promoter, such as heat shock promoter, may be used in this invention.

The super-compound interferon may also be produced with its gene as artificially synthesized cDNA with adjustment of its sequence from the wild-type according to codon preference of E. Coli. Extensive discussion of said codon usage (preference) may be found in U.S. Patent No. 4,695,623. See e.g. column 6, line 41 - column 7, line 35

The above described super-compound interferon possesses anti-viral or anti-tumor activity and therefore useful in preventing and treating viral diseases, tumors or cancers.

The virus diseases include but are not limited to Hepatitis A, Hepatitis B, Hepatitis C, other types of hepatitis, infections caused by Epstein-Barr virus, Cytomegalovirus, herpes simplex viruses, other herpes viruses, papovaviruses, poxviruses, picornaviruses, adenoviruses, rhinoviruses, human T cell leukaemia viruses I, human T cell leukaemia viruses II, or human T cell leukemia viruses III.

Therefore, this invention provides a method for inhibitive virus replication or virus infected cells by contacting said virus or infected cells with an effective amount of the super-compound interferon or its equivalent. This super-compound interferon is useful in preventing or treating the following cancers or tumors:

Cancer	Skin Cancer	Basal Cell Carcinomas
		Malignant Melanoma

	Renal cell carcinoma	
	Liver Cancer	
	Thyroid Cancer	
	Rhinopharyngeal Cancer	
	Solid Carcinoma	Prostate Cancer
		Tummy Cancer
		Esophagus Cancer
		Recta Cancer
		Pancreas Cancer
		Mammary Cancer
	Ovarian Cancer & Superficial Bladder Cancer	
	Hemangioma	
	Epidermoid Carcinoma	Cervical Cancer
		Not-small Cell Lung Cancer
		Small Cell Lung Cancer
		Glioma
Malignant Hemal Disease	Leucocythemia	Acute Leucocythemia
		Chronic Leucocythemia
	Chronic Myelocytic Leukemia	
	Hairy Cell Leukemia	
	Lymphadenoma	
	Multiple Myeloma	
Others	Polycythemia Vera	
	Kaposi's Sarcoma	

Accordingly, this invention provides a method for inhibiting tumor or cancer cell growth by contacting the super-compound interferon or its equivalent with said tumor or cancer cells. In a further embodiment, the super-compound interferon inhibits the DNA duplication and secretion of HBsAg and HBeAg of Hepatitis B Virus.

This invention also provides an artificial gene codes for the super-compound interferon or its equivalent. It is within the ordinary skill to design an artificial gene. Many methods for generating nucleotide sequences and other molecular biology techniques have been described previously. See for example, Joseph Sambrook and David W. Russell, Molecular Cloning: A laboratory Manual, December 2000, published by Cold Spring Harbor Laboratory Press.

This invention provides a vector comprising the gene which codes for the super-compound interferon or its equivalent.

This invention provides an expression system comprising the vector comprising the gene which codes for the super-compound interferon or its equivalent. The cells include but are not limited to prokaryotic or eukaryotic cells.

This invention also provides a host cell comprising the vector comprising the gene which codes for the super-compound interferon or its equivalent.

This invention provides a process for production of recombinant super-compound interferon comprising introducing an artificial gene with selected codon preference into an appropriate host, culturing said introduced host in an appropriate condition for the expression of said compound interferon and harvesting the expressed compound interferon.

The process may comprise extraction of super-compound interferon from fermentation broth, collection of inclusion body, denaturation and renaturation of the harvested protein.

The process may maintain the high efficacy even when the super-compound interferon is used with an agent and in a particular concentration. The process also comprises separation and purification of the super-compound interferon. The process further comprises lyophilization of the purified super-compound interferon. The process comprises production of liquid injection of super-compound interferon.

This invention also provides the produced super-compound interferon by the above processes.

This invention provides a composition comprising the recombinant super-compound interferon or its equivalent and a suitable carrier.

This invention provides a pharmaceutical composition comprising the recombinant super-compound interferon or its equivalent and a pharmaceutically acceptable carrier.

This invention provides a method for treating viral diseases or tumor in a subject comprising administering to the subject an effective amount of the super-compound interferon or its equivalent.

This invention provides the above-described method wherein the viral diseases is Hepatitis A, Hepatitis B, Hepatitis C, other types of hepatitis, infections of viruses caused by Epstein-Barr virus, Cytomegalovirus, herpes simplex viruses, or other type of herpes viruses, papovaviruses,

poxviruses, picornaviruses, adenoviruses, rhinoviruses, human T cell leukaemia viruses I, or human T cell leukaemia viruses II, or human T cell leukemia virus III.

5 This invention provides the above-described method wherein super-compound interferon was administered via oral, vein injection, muscle injection, peritoneal injection, subcutaneous injection, nasal, mucosal administration, by
10 inhalation via an inspirator.

This invention provides the above-described method wherein super-compound interferon was administered following the protocol of injection 9 μ g or 15 μ g per day, 3 times a
15 week, totally 24 weeks.

It was surprising to find that rSIFN-co, the spatial configuration of which has been changed, is not only a preparation to inhibit the DNA duplication of Hepatitis B,
20 but to inhibit the secretion of HBsAg and HBeAg.

One objective of this invention is to offer a preparation of rSIFN-co to directly inhibit the DNA duplication of Hepatitis B viruses and the secretion of HBeAg and HBsAg of
25 Hepatitis B and decrease them to normal levels.

In one of the results of this invention, rSIFN-co was produced with recombinant techniques. On the condition of fixed amino acid sequence, the IFN DNA was redesigned
30 according to the *E. Coli*. codon usage and then the rSIFN-co gene was artificially synthesized. rSIFN-co cDNA was cloned into the high-expression vector of *E. Coli*. by DNA recombinant techniques, and a high expression of rSIFN-co was gained by using of induce/activate-mechanism of L-arabinose to activate the transcription of P_{BAD} promoter.
35

Compared with usual thermo-induction, pH induction and IPTG

induction systems of genetic engineering, arabinose induction/activation system has some advantages: (1) Common systems relieve promoter function by creating a "derepression" pattern. Promoters then induce downstream gene expression. So temperature and pH change and the addition of IPTG cannot activate promoters directly. In the system disclosed herein, L-arabinose not only deactivates and represses but also activates the transcription of *P_{lac}* promoter which induces a high expression of rSIFN-co. Therefore, the arabinose induction/activation system is a more effective expression system. (2) The relation between Exogenous and L-arabinose dosage is linearity. This means the concentration of arabinose can be changed to adjust the expression level of the exogenous gene. Therefore, it is easier to control the exogenous gene expression level in *E. Coli.* by arabinose than by changing temperature and pH value. This characteristic is significant for the formation of inclusion bodies. (3) L-arabinose is resourceful cheap and safe, which, on the contrary, are the disadvantages of other inducers such as IPTG.

This embodiment creates an effective and resistant rSIFN-co-expressing *E. Coli.* engineering strain with an L-arabinose induction/activation system. The strain is cultivated and fermented under suitable conditions to harvest the bacterial bodies. Inclusion bodies are then purified after destroying bacteria and washing repeatedly. The end result, mass of high-purity, spacial-structure-changed rSIFN-co protein for this invention and for clinical treatment, was gained from denaturation and renaturation of inclusion bodies and a series of purification steps.

The following are some rSIFN-co preparations: tablets, capsules, oral liquids, pastes, injections, sprays, suppositories, and solutions. Injections are recommended.

It is common to subcutaneously inject or vein-inject the medicine. The medicine carrier could be any acceptance medicine carrier, including carbohydrate, cellulose, adhesives, disintegration agents, emollients, filling, add-
5 dissolve agent, amortization, preservative, add-thick agent, matching, etc.

This invention also provides a pharmaceutical composition comprising the above composition and a pharmaceutically
10 acceptable carrier.

For the purposes of this invention, "pharmaceutically acceptable carriers" means any of the standard pharmaceutical carriers. Examples of suitable carriers are
15 well known in the art and may include, but are not limited to, any of the standard pharmaceutical carriers such as a phosphate buffered saline solution and various wetting agents. Other carriers may include additives used in
tablets, granules and capsules, etc. Typically such
20 carriers contain excipients such as starch, milk, sugar, certain types of clay, gelatin, stearic acid or salts thereof, magnesium or calcium stearate, talc, vegetable fats or oils, gum, glycols or other known excipients. Such carriers may also include flavor and color additives or
25 other ingredients. Compositions comprising such carriers are formulated by well-known conventional methods.

This invention will be better understood from the examples which follow. However, one skilled in the art will readily
30 appreciate that the specific methods and results discussed are merely illustrative of the invention as described more fully in the claims which follow thereafter.

EXPERIMENTAL DETAILS

EXAMPLE 1

rSIFN-co is a new interferon molecule constructed according to conservative amino acid in human IFN- α subtype with genetic engineering method. It has been proven that rSIFN-co has broad-spectrum IFN activity, such as high antiviral and tumor inhibition activity, especially for effectively treating Hepatitis C.

E. Coli. codon was used to redesign rSIFN-co cDNA and then artificially synthesize cDNA of rSIFN-co from published rSIFN-co DNA sequences and deduced amino acid sequences (Figure 1).

In order to get pure rSIFN-co protein, rSIFN-co cDNA was cloned into E. Coli. high-expression vector, and L-arabinose, which can activate strong PBAD promoter in vectors, was used to induce high expression of rSIFN-co gene.

1. Synthesis of E. Coli. cDNA Sequence

1.1 Redesign of rSIFN-co cDNA sequence

rSIFN-co cDNA was redesigned according to the codon usage of E. Coli. to achieve high expression in E. Coli. Deduced amino acid sequence from the redesigned cDNA sequence of rSIFN-co is completely coincidental with primitive amino acid sequence of published rSIFN-co (Figure 1).

1.2 rSIFN-co cDNA sequence synthesis

1.2.1 rSIFN-co cDNA 5'-terminus and 3'- terminus semi-molecular synthesis

Two semi-moleculars can be directly synthesized: rSIFN-co cDNA 5'- terminus 280bp (fragment I) and 3'- terminus 268bp(fragment II) by PCR. There are 41bp overlapping among fragment II and fragment I.

(1) Chemical synthesis oligodeoxynucleotide fragment:

Oligomer A:

5'ATCTGGGACCTGCCGACAGCCCACTCCCTGGGTACCGTCTGCTCTGATCCTGCTGCTCA
GATGCGTCTGATCTCTCCCGTTCTCTGCTGAAAGACCGTCAAGAC3'

Oligomer B:

5'CTGAAAGACCGTCACGACTTCGGTTCCCGCAGGAGAGGTTCCAGCGTAACAGTTCCAGA
AGCTCAGGCTATCTCCGTTCTGCACGAAATGATCCAGCAGACCTTC3'

Oligomer C:

5 GCTGCTGGTACAGTTCGGTGTAGAATTTTCCAGCAGGAGTTGGTCCCAAGCAGCGGAGGAG
TCATTGGTGGAGACAGGTTGAGGCTCTGCTGGATCATTTTC3'

Oligomer D:

5'ATCCCTGCTGGAAAAATTTCTACACCGAACTGTACACAGCAGCTGAACGAACTGGAAGCTTGGG
TTATCCAGGAAGTTGGTGTGAGAAACCCCGCTGATGAAC3'

Oligomer E:

5'GAAGAAACCCCGCTGATGAACGTTGACTCCATCCTGGCTGTAAAAAATACTTCCAGCGTAT
CACCTGTACTTGACCGAAAAAATACTCCCGTGGCTTGGG3'

Oligomer F:

5'TTATTCCTTACGACGACAGCTTCCTGCAGGTTGGTGGACAGGGAGAAGGAACGCATGATT
13 CAGCAGAACAACTTCCCAAGGGCACGGGAGATTTTTTTTCGGTCAGG3'

PCR I for Fragment I: oligodeoxynucleotide B as template,
oligodeoxynucleotide A and C as primers, synthesized 280 bp
Fragment I.

20 PCR I mixture (units: μ l)

sterilized distilled water	39
10 \times Pfu buffer (Stratagen American Ltd.)	5
dNTP mixture (dNTP concentration 2.5 mmol/L)	2
Oligomer A primer (25 μ mol/L)	1
Oligomer C primer (25 μ mol/L)	1
Oligomer B template (1 μ mol/L)	1
Pfu DNA polymerase (Stratagen American Ltd.) (25 U/ μ l)	1
Total volume	50 μ l

PCR cycle: 95 I

2m-(95°C45s-65°C1m-72°C1m) \times 25 cycle-72°C10m-4°C

25

PCR II for Fragment II: oligodeoxynucleotide E as
template, oligodeoxynucleotide D and F as primers,
synthesized 268bp Fragment II.

PCR II mixture (units: μ l)

sterilized distilled water	39
10xPfu buffer (Stratagen American Ltd.)	5
dNTP mixture (dNTP concentration 2.5mmol/L)	2
Oligomer D primer (25 µmol/L)	1
Oligomer F primer (25 µmol/L)	1
Oligomer E template (1 µmol/L)	1
Pfu DNA polymerase (Stratagen American Ltd.) (25U/µl)	1
Total volume	50µl

PCR cycle: the same as PCR I

1.2.2 Assembling of rSIFN-co cDNA

- 5 Fragment I and II were assembled together to get the complete cDNA molecular sequence of rSIFN-co using the overlapping and extending PCR method. Restriction enzyme Nde I and Pst I were introduced to clone rSIFN-co cDNA sequence into plasmid.

10 (1) Chemical synthesis primers

Oligomer G: 5'ATCGCCATATGTGCGACCTGCGCAGACCC3'

Oligomer H: 5'ACTGCCAGGCTGCAGTTATCTTTACGACGACAGGTTCC3'

(2) Overlapping and extending PCR

	(units: µl)
sterilized distilled water	38
10xPfu buffer (Stratagen American Ltd.)	5
dNTP mixture (dNTP concentration 2.5mmol/L)	2
primer G (25 µmol/L)	1
primer H (25 µmol/L)	1
*fragment I production (1 µmol/L)	1
*fragment II production (1 µmol/L)	1
Pfu DNA polymerase (Stratagen American Ltd.) (2.5U/µl)	1
15 Total volume	50µl

*Separate and purify PCR production with StrataPrep PCR purification kit produced by Stratagen American Ltd. And dissolve into sterilized distilled water.

PCR cycle: the same as PCR I

2. rSIFN-co gene clone and sequence analysis

pLac T7 plasmid as cloning vector. pLac T7 plasmid is

- 5 reconstructed with pBluescript II KS(+) plasmid produced by Stratagen (Figure 3).

Purified PCR production of rSIFN-co cDNA with StrataPrep

PCR purification kit. Digest cDNA and pLac T7 plasmid with

- 10 NdeI and PstI. Run 1% agarose gel electrophoresis and separate these double-digested DNA fragments. Recover 507bp long rSIFN-co DNA fragment and 2.9kb plasmid DNA fragment.

Ligate these fragments by T4 DNA ligase to form a recombinant plasmid. Transform DH₅ α competent cells (Gibco) with the recombinant plasmid. culture at 37°C overnight.

- 15 Identify the positive recombinant colony, named PHY-1.

Run DNA sequencing with SequiTherm™ Cycle Sequencing Kit

produced by American Epicentre Technologies Ltd using Li-

COR Model 4000L. Primers are T7 and T3 common sequence

- 20 primer, the DNA sequencing result matches theoretic design. Purify the rSIFN-co, sequence the N-terminus amino acids, the N-terminus amino acid sequence matches experimental design which is as follows:

N- Cys-Asp-Leu-Pro-Gln-Thr-His-Ser-Leu-Gly-Asn-Arg-Arg-Ala-Leu-

- 25 3. Construction, transformation, identification, and hereditary stability of expression vector

3.1 Construction and transformation of expression vector

Digested E. Coli. expression vector PHY-4 (see Figure 3)

with Nde I to linearize and subsequently digest with Xba I.

- 30 Run 1% agarose gel electrophoresis, and purify the 4.8kb PHY-4 Nde I-Xba I digest fragment with QIAEX II kit produced

18

by QIAGEN Germany Ltd.

At the same time, the pHY-4 plasmid is double digested with Nde I-Xba I. Run 1% agarose gel electrophoresis and purify the 715bp fragment. Ligate the rSIFN-co and pHY-4 fragments with T4 DNA ligase to construct the recombinant plasmid (See Figure 4). Transform DH₅ competent cells with the recombinant plasmid. Spread the transformed cells on LB plate with Amp, 37°C culture overnight.

3.2 Positive cloning strain screening

- Randomly choose *E. Coli.* colonies from above LB-plate, screening the positive strains containing recombinant vector by endonuclease digesting and PCR analysis. Name one of the positive recombinant plasmid pHY-5, and name the strain containing pHY-5 plasmid PVIII. Amplify and store the positive strain with glycerol in -80°C.

4. High expression of rSIFN-co gene in *E. Coli.*

- In pHY-5 plasmid, rSIFN-co gene is under control of strong promoter P_{BAD}. This promoter is positively and negatively regulated by the product of the gene *araC*. AraC is a transcriptional regulator that forms a complex with arabinose. In the absence of arabinose, the AraC dimer binds O₂ and I₁ forming a 210bp loop. This conformation leads to a complete inhibition of transcription. In the presence of arabinose, the dimer is released from O₂ and binds I₁ and I₂ leading to transcription. Arabinose binding deactivates, represses and even activates the transcription of P_{BAD} promoter, which stimulates P_{BAD} inducing high expression of rSIFN-co. rSIFN-co expression level in PVIII is more than 50% of the total *E. Coli.* proteins.

5. Summary

rSIFN-CO is a new interferon molecule artificially built

according to the conservative amino acid of human α interferons. It has been proven as a effective anti-hepatitis drug. In order to get enough pure rSIFN-co protein, a stable recombinant *E. Coli.* strain which high expresses rSIFN-co protein was constructed.

First, according to published rSIFN-co amino acid sequence, *E. Coli.* codon was used to synthesize whole cDNA of rSIFN-co. This DNA fragment was sequenced and proved that the 501bp codon sequence and TAA termination codon sequence are valid and identical with theocratic design. Subsequent analysis revealed that the N-terminus amino acid sequence and amino acid composed of rSIFN-co produced by the recombinant strain were both identical to the prediction.

The rSIFN-co cDNA was cloned into *E. Coli.* high-expression vector pHY-4 plasmid to construct the recombinant plasmid pHY-5. *E. Coli.* LMG194 strain was further transformed with pHY-4 plasmid to get stable rSIFN-co high-expression transformant. This transformant was cultured for 30 generations. The heredity of pHY-5 recombinant plasmid in *E. Coli.* LMG194 was normal and stable, and the expression of rSIFN-co was high and steady.

E. Coli. LMG194, which contains recombinant pHY-5 plasmid, is actually an ideal high-expression engineering strain.

6. References

- 1). Blatt LM, Davis JM, Klein SB. et al. The biologic activity and molecular characterization of a novel synthetic interferon-alpha species, consensus interferon. *Journal of Interferon and Cytokine Research*, 1996;16(7):489-499.
- 2). Alton, K. et al: Production characterization and biological effects of recombinant DNA derived human IFN- α

- and IFN- γ analogs. In: De Maeger E, Schellekens H. eds. The Biology of Interferon System. 2nd ed. Amsterdam: Elsevier Science Publishers, 1983: 119-128
- 3). Pfeffer LM. Biologic activity of natural and synthetic type 1 interferons. *Seminars in Oncology*, 1997;24 (3 suppl 9):S9-63--S9-69.
- 4). Ozes ON, Reiter Z, Klein S, et al. A comparison of interferon- α with natural recombinant interferons- α : antiviral, antiproliferative, and natural killer-inducing activities. *J. Interferon Res.*, 1992; 12:55-59.
- 5). Heathcote EJJ, Keefe EB, Lee SS, et al. Re-treatment of chronic Hepatitis C with consensus interferon. *Hepatology*, 1998;27(4):1136-1143.
- 6). Klein ML, Bartley TD, Lai PH, et al. Structural characterization of recombinant consensus interferon- α . *Journal of Chromatography*, 1988; 454:203-215.
- 7). The Wisconsin Package, by Genetics Computer Group, Inc. Copyright 1992, Madison, Wisconsin, USA
- 8). Nishimura, A et al: A rapid and highly efficient method for preparation of competent *E. coli* cells. *Nucleic Acids Res.* 1990, 18:6169
- 9). All molecular cloning techniques used are from Sambrook, J., E. F. Fritsch and T. Maniatis. *Molecular Cloning: A laboratory manual* 2nd ed. CSH Laboratory Press, Cold Spring Harbour, NY. 1989.
- 10). Guzman, L. M et al: Tight regulation, modulation, and high-level expression by vectors containing the arabinose *P_{BAD}* promoter. *J. Bacteriol.* 1995, 177 4121- 4130.

30 EXAMPLE 2

Separation and purification of rSIFN- α

1. Fermentation

Inoculate the recombinant strain in LB media, shaking (200 rpm) under 37°C overnight (approximate 18 h), then add 30%

glycerol to the fermentation broth to get final concentration of 15%, allotted to 1 ml tube and kept in -20°C as seed for production.

- 5 Add 1% of the seed to LB media, shaking (200 rpm) under 37°C overnight to enlarge the scale of the seed, then add to RM media with a ratio of 10%, culturing under 37°C. Add arabinose (20% solution) to 0.02% as an inductor when the OD600 reaches about 2.0. 4 hours after that, stop the culture process, collect the bacteria by centrifuge, resuspend the pellet with buffer A, and keep in -20°C overnight. Thaw and break the bacteria by homogenizer, then centrifuge. Wash the pellet with buffer B, buffer C, and distilled water to get a relatively pure inclusion body.

15

2. Denaturation and renaturation

- Dissolve the inclusion body in Guanidine-HCl (or urea) of 6 mol/L. The solution will be a little cloudy. Centrifuge it at a speed of 10000 rpm. Determines the protein concentration of the supernatant. This supernatant is called "denaturation solution." Add the denaturation solution to renaturation buffer, and keep the final protein concentration under 0.3 mg/ml. It is better to add the totally denaturation solution in three steps instead of one step. Keep the solution overnight under 4°C. Afterwards, dialyze against 10 mol/L and 5 mol/L PB buffer and distilled water, then adjusting its pH by 2 mol/L HAc-NaAc. Let it standstill for a while, then filtrate.

30

3. Purification

POROS HS/M anion exchange chromatography:

Equivalent column with 20 mmol/L HAc-NaAc(pH 5.0)



35

Load samples at a speed of 30 ml/min

↓
Wash with 20 CV 20 mmol/L HAc-NaAc(pH 5.0)
↓
5 5 CV of 0.15 mol/L NaCl+20 mmol/L HAc-NaAc(pH 5.0) wash
↓
3 CV of 0.18 mol/L NaCl+20 mmol/L HAc-NaAc(pH 5.0) wash
10 ↓
0.25 mol/L NaCl + 20 mmol/L HAc-NaAc(pH 5.0) elute target
protein
↓
15 Chelating sepharose™ fast flow:
Add PB buffer of 0.2 mol/L (pH 6.6) and NaCl of 4 mol/L in
the solution from HS to adjust solution pH to pH 6.0 and
NaCl concentration to 1 mol/L.
equivalent Column with buffer D
20 ↓
Loading at a rate of 1 ml/min
↓
25 Wash with buffer E
↓
Wash with buffer F
↓
30 Elute with buffer G
↓
Condense the eluted solution by POROS HS/M. Sometimes a
step of purification by sephacryl S-100 can be added to
35 meet with stricter purity requirements.

↓
Wash with 20 CV 20 mmol/L HAc-NaAc (pH 5.0)
↓
5 CV of 0.15 mol/L NaCl + 20 mmol/L HAc-NaAc (pH 5.0) wash
↓
3 CV of 0.18 mol/L NaCl + 20 mmol/L HAc-NaAc (pH 5.0) wash
↓
0.25 mol/L NaCl + 20 mmol/L HAc-NaAc (pH 5.0) elute target
protein
↓
Chelating sepharose™ fast flow:
Add PB buffer of 0.2 mol/L (pH 6.6) and NaCl of 4 mol/L in
the solution from HS to adjust solution pH to pH 6.0 and
NaCl concentration to 1 mol/L.
equivalent Column with buffer D
↓
Loading at a rate of 1 ml/min
↓
Wash with buffer E
↓
Wash with buffer F
↓
Elute with buffer G
↓
Condense the eluted solution by POROS HS/M. Sometimes a
step of purification by sephacryl S-100 can be added to
meet with stricter purity requirements.

Note:

- Buffer A: 100 mmol/L Tris-HCl, pH 7.5-10 mmol/L EDTA-100 mmol/L NaCl
- 5 Buffer B: 50 mmol/L Tris-HCl, pH 7.5-1 mol/L Urea-10 mmol/L EDTA-0.5% Triton X-100
- Buffer C: 50 mmol/L Tris-HCl, pH 7.5-2 mol/L Urea-10 mmol/L EDTA-0.5% Triton X-100
- 10 Buffer D: 1 mol/L NaCl ---50 mmol/L Na_2HPO_4 (pH 5.5)
- Buffer E: 1 mol/L NaCl ---50 mmol/L Na_2HPO_4 (pH 5.0)
- Buffer F: 1 mol/L NaCl ---50 mmol/L Na_2HPO_4 (pH 4.0)
- Buffer G: 1 mol/L NaCl ---50 mmol/L Na_2HPO_4 (pH 3.6)
- Renaturation buffer: 0.5 mol/L Arginine-150 mmol/L Tris-HCl, pH 7.5-0.2 mmol/L EDTA
- 15 LB Media: 1 L
- Tryptone 10 g
- Yeast extracts 5 g
- NaCl 10 g
- RM Media: 1 L
- 20 Casein 20 g
- MgCl 1 mmol/L (0.203 g)
- Na_2HPO_4 4 g
- KH_2PO_4 3 g
- NaCl 0.5 g
- 25 NH_4Cl 1 g

After purification, the buffer was changed to PBS (pH 7.0) along with the step of condensing by POROS HS/M. This is called the "Protein Stock Solution." It can be directly used in the preparation of injections or sprays, or stored at 2-8 degree centigrade.

30

Formula for injection:

	Solution	Lyophilized powder
Solution of rSIFN-co	34.5 µg/ml	34.5 µg/ml
FB (pH7.0)	25mmol/L	10mmol/L
Glycine	-----	0.4mol/L
NaCl	0.1mol/L	-----

For spray:

EDTA	0.01%
Tween 80	0.05%
Trisodium citrate	10mmol/L
Glycerol	1.26%
Sodium Chloride	0.03%
Phenylmethanol	0.5%
HSA	0.1%
rSIFN-co	10 µg/ml

5 QUALITY CONTROL PROCESS

During purification tests for protein content, protein purity, specific activity and pyrogen are conducted after each step. When the stock solution is obtained, all the tests listed in the table are done one after the other.

10

The quality of the product is controlled according to "Chinese Requirements for Biologics"

1. Original protein solution

Lowry

Item of Test	Method
Protein Stock Solution:	
Test for Protein Content	Lowry
Test for Protein Purity	Non-reductive SDS-PAGE (sodium dodecyl sulfate polyacrylamide gel electrophoresis) HPLC Analysis
Test for Molecular Weights	Reductive SDS-PAGE

Test for Specific Activity	According to Method in "Specific Activity Test of Interferon"
Test for Leftover Exogenous DNA	Using DNA Labeling and Detection Kit
Test for Activity of Leftover Antibiotics	According to Method in "Chemical and Other Test Methods for Biologics"
Test for Bacterial Endotoxin	According to Method in "Requirements for Bacterial Endotoxin Test of Biologics"
Test for Isoelectric Point	Isoelectric Focusing Electrophoresis
Test for Identify Characteristics of the Protein	UV spectrum (range of wavelength: 190-380nm)
	Peptide Mapping (hydrolyzed by pancreatic enzyme, analyzed by C-18 column)
	N-terminal Sequence Test
	C-terminal Sequence Test
	Circular Dichroism
	Amino Acid Analysis
Semi-finished Product	
Test for Bacterial Endotoxin	According to Method in "Requirements for Bacterial Endotoxin Test of Biologics"
Product	
Appearance Check	
Chemical	According to Method in "Chemical and Other Test Methods for Biologics"
Test for Specific Activity	According to Method in "Specific Activity Test of Interferon"
Sterility Test	According to Method in "c"
Abnormal Toxicity Test	Test on Mouse
Pyrogen Test	According to Method in "Requirements for Pyrogen Test of Biologics"
Test for Stability of Product	

Note: "Chemical and Other Test Methods for Biologics", "Requirements for Pyrogen Test of Biologics" and "Requirements for Bacterial Endotoxin Test of Biologics" all can be found in the "Chinese Requirements for Biologics." "Chinese Requirements for Biologics," PAN

Zhengan, ZHANG Xinhui, DUAN Zhibing, et al. Chinese Biologics Standardization committee. Published by Chemical Industry Publishing Company, 2000.

5 **EXAMPLE 3**

Stability of lyophilized Powder of Recombinant Super-Compound Interferon Injection

10 The stability experiments were carried out with samples of lyophilized powder of recombinant super-compound interferon (rSIFN-co) injection in two specifications and three batches. The experiments started on April, 2000.

1. Sample Source

15 Samples were supplied by Sichuan Huiyang Life-engineering Ltd., Sichuan Province. Lot: 990101-03, 990101-05, 990102-03, 990102-05, 990103-03, 990103-05

2. Sample Specifications

20 Every sample in this experiment should conform with the requirements in the table below.

Table 1 Standards of Samples in Experiment

Items	Standards
1. Appearance	white loose powder
2. Dissolving time	dissolve rapidly in injection water(within 2 min) at room temperature
3. Clarity	colourless liquid or with little milk-like glisten; should not be cloudy, impurity or with indiscerptible deposit
4. pH value	6.5-7.5
5. Potency (IU/dose)	80%-150% of indicated quantity (9 μ g:4.5 * 10 ⁶ IU, 15 μ g: 7.5 * 10 ⁶ IU)
6. Moisture	no more than 3.0% (W/W)

3. Experiment Content

15.3.1 Test samples at 2-8°C: The test samples were put into a 2-8°C refrigerator, then the above items of these samples were respectively tested in the 1st, 3rd, 6th, 9th, 12th, 18th, 24th, 30th, 36th month. The results were recorded.

15.3.2 Test samples at 25°C: The test samples were put into a thermostat at 25°C, then the above items of these samples were respectively tested in the 1st, 3rd, 6th, 9th, 12th, 18th, 24th, 30th month. The results were recorded.

15.3.3 Test samples at 37°C: The test samples were put into a thermostat at 37°C, then the above items of these samples were respectively tested in the 1st, 3rd, 6th, 9th, 12th, 18th, 24th month. The results were recorded.

4. Results and Conclusion

1) At 37°C, according to data collected at designated points during testing and compared with data before testing, the potency began descending from the 6th month and the changes in the three batches were similar. The appearance of other items had no changes.

2) At 25°C, according to data collected at designated points during testing and compared with data before the testing, the potency only had a little change, and the changes in the three batches were similar. The appearance of other items had no changes.

3) At 2-8°C, according to data collected at designated points during testing and compared with data before testing, the potency of the three batches all were stable. The appearance of other items also had no changes.

In conclusion, it is suggested that the lyophilized powder of recombinant super-compound interferon for injection should be better stored and transported at low temperatures. Without such conditions, the product can also be stored for short periods (i.e. 3 months) at room

temperature.

EXAMPLE 4

5 rSIFN- α 2b inhibits HBV-DNA duplication and secretion of
HBsAg and HBeAg.

Materials

10 Solvent and Dispensing Method: Add 1ml saline into each
vial, dissolve, and mix with MEM culture medium at
different concentrations. Mix on the spot.

Control drugs: IFN- α 2b (Intron A) as lyophilized powder,
purchased from Schering Plough. 3×10^6 U each, mix to
15 3×10^6 U/ml with culture medium; INFERGEN[®] (Liquid solution)
, purchased from Amgen, 9 μ g, 0.3ml each, equal to 9×10^6 U,
and mix with 9×10^6 U/ml culture medium preserves at 4°C;
2.2.15 cell: 2.2.15 cell line of hepatoma (Hep G2) cloned
and transfected by HBV DNA, constructed by Mount Sinai
Medical Center.

20 Reagent: MEM powder, Gibco American Ltd. cattle fetal blood
serum, HycloneLab American Ltd. G-418(Geneticin); MEM
dispensing, Gibco American Ltd.; L-Glutamyl, imported and
packaged by JING KE Chemical Ltd.; HBsAg and HBeAg solid-
25 phase radioimmunoassay box, Northward Reagent Institute of
Chinese Isotope Ltd.; Biograngetina, Northern China
Medicine; And Lipofectin, Gibco American Ltd.

Experimental goods and equipment: culture bottle, Denmark
30 TuncionTM, 24-well and 96-well culture board, Corning
American Ltd.; Carbon Dioxide hatching box, Shel-Lab
American Ltd.; MEM culture medium 100ml; 10% cattle fetal
blood serum, 3% Glutamyli¹, G418 380 μ g/ml,
biograngetina 500U/ml.

35 Method:

2.2.15 cell culture: Added 0.25% pancreatic enzymes into culture box with full of 2.2.15 cell, digest at 37°C for 3 minutes, and add culture medium to stop digest and disturb it to disperse the cells, reproduce with ratio of 1:3. They will reach full growth in 10 days.

Toxicity test: Set groups of different concentrations and a control group in which cell is not acted on with medicine. Digest cell, and dispense to a 100,000 cell/ml solution. Inoculate to 96-well culture board, 200µl each well, culture at 37°C for 24h with 5% CO₂. Test when simple cell layer grows.

Dispense rSIFN-co to 1.8×10⁷ IU/ml solution than prepare a series of solutions diluted at two-fold gradients. Add into 96-well culture board, 3 wells per concentration. Change the solution every 4 days. Test cytopathic effect by microscope after 8 days. Fully destroy as 4, 75% as 3, 50% as 2, 25% as 1, zero as 0. Calculate average cell lesion and inhibition rate of different concentrations. Calculate TC50 and TCO according to the Reed Muench method.

$$TC50 = \text{Antilog} \left(B + \frac{50-B}{A-B} \times C \right)$$

A=log >50% medicine concentration, B=log<50% medicine concentration, C=log dilution power

Inhibition test for HBsAg and HBsAg: Separate into positive and negative HBsAg and HBsAg contrast groups, cell contrast group and medicine concentration groups. Inoculate 700,000 cells/ml of 2.2.15 cell into 6-well culture board. 3 ml each well, culture at 37°C for 24h with 5% CO₂, then prepare 5 gradiently diluted solutions with 3-fold as the grade (Prepare 5 solutions, each with a different protein concentration. The concentration of Solution 2 is 3 times lower than that of Solution 1, the concentration of Solution 3 is 3 times lower than that of Solution 2, etc.)

4.5*10⁶IU/ml, 1.5*10⁶IU/ml, 0.5*10⁶IU/ml, 0.17*10⁶IU/ml, and 0.056*10⁶IU/ml, 1 well per concentration, culture at 37°C for 24h with 5% CO₂. Change solutions every 4 days using the same solution. Collect all culture medium on the 6th day.

- 5 Preserve at -20°C Repeat test 3 times to estimate HBsAg and HBeAg with solid-phase radioimmunoassay box (Northward Reagent Institute of Chinese Isotope Ltd.). Estimate cpm value of each well with a γ- accounting machine.

- 10 Effects calculation: Calculate cpm mean value of contrast groups and different-concentration groups and their standard deviation, P/N value such as inhibition rate, IC50 and SI.

$$1) \text{ Antigen inhibition rate (\%)} = \frac{A-B}{A} \times 100$$

- 15 A = cpm of control group; B = cpm of test group;

- 2) Counting the half-efficiency concentration of the medicine

$$\text{Antigen inhibition IC50} = \text{Antilog} \left(B + \frac{50-B}{A-B} \times C \right)$$

- 20 A=log>50% medicine concentration, B=log<50%medicine concentration, C=log dilution power

- 3) SI of interspace-conformation changed rSIFN-co effect on HBeAg and HBsAg in 2.2.15 cell culture:

$$SI = \frac{IC50}{IC50}$$

- 25 4) Estimate the differences in cpm of each dilution degree from the control group using student t test

- 30 Southern blot: (1) HBV-DNA extract in 2.2.15 cell: Culture cell 8 days. Exsuction culture medium (Separate cells from culture medium by means of draining the culture medium.). Add lysis buffer to break cells, then extract 2 times with a mixture of phenol, chloroform and isoamyl alcohol

(1:1:1), 10,000g centrifuge. Collect the supernatant adding anhydrous alcohol to deposit nucleic acid. Vacuum draw, redissolve into 20 μ l TE buffer. (2) Electrophoresis: Add 6X DNA loading buffer, electrophoresis on 1.5% agarose gel, IV/cm, at fixed pressure for 14-18h. (3) Denaturation and hybridization: respectively dip gel into HCl, denaturation buffer and neutralization buffer. (4) Transmembrane: Make an orderly transfer of DNA to Hybond-N membrane. Bake, hybridize and expose with dot blot hybridization. Scan and analyze relative density with gel-pro software. Calculate inhibition rate and IC50.

Results

Results from Tables 1, 2 and 3 show: After maximum innocuous concentration exponent culturing for 8 days with 2.2.15 cell, the maxima is $9.0 \pm 0 \times 10^4$ IU/ml average inhibition rate of maximum innocuous concentration rSIFN-co to HBeAg is $46.0 \pm 5.25\%$ ($P < 0.001$), IC50 is $4.54 \pm 1.32 \times 10^4$ IU/ml, SI is 3.96; rate to HBeAg is $44.8 \pm 6.6\%$, IC50 is $6.49 \pm 0.42 \times 10^4$ IU/ml, SI is 2.77. This shows that rSIFN-co can significantly inhibit the activity of HBeAg and HBeAg; but that the IFN of the contrast group and INFERGEN[®] cannot. It has also been proved in clinic that rSIFN-co can decrease HBeAg and HBeAg or return them to normal levels.

Table 1: Results of inhibition rate of rSTPn-co to HBsAg and HBsAg
First batch: (rSTPn-co)

				Inhibition effect to HBsAg				Average inhibition rate	Accumulation n	1- n	Accumulated inhibition rate
Concentratio n(×10 ⁴ U/ml)	First well	Second well	Third well	First well	Second well	Third well	Fourth well				
900	9026	8976	8947	0.436227	0.439315	0.44565	0.407079	0.44565	0.392921	0.392921	0.614653416
300	9616	12082	1008	0.399375	0.24534	0.36926	0.37987	0.36926	0.398289	1.254924	0.301392321
100	9822	16002	1280	0.36508	0.0005	0.2005	0.195826	0.2005	0.300833	2.95008	0.09867188
33.3333	1577	19306	0	0.014991	0	0	0.064997	0.064997	0.0049569	3.451091	0.001633453
11.1111	1917	22230	1893	0	0	0	0	0	0	4.054091	0
Control	Cell 16010			Blank	0	0	0	0	0	IC50	602.7444601
				Inhibition effect to HBsAg				Average inhibition rate	Accumulation n	1- n	Accumulated inhibition rate
Concentratio n(×10 ⁴ U/ml)	First well	Second well	Third well	First well	Second well	Third well	Fourth well				
900	7706	7240	7134	0.342155	0.38193	0.39269	0.37261	0.392258	0.627759	0.392258	0.595006426
300	8856	7770	9476	0.243901	0.33600	0.19305	0.257014	0.5499972	1.370724	0.286319225	0.286319225
100	1081	10720	1033	0.07649	0.03485	0.11814	0.093165	0.292883	2.27756	0.113972619	0.113972619
33.3333	1074	11114	1057	0.082807	0.05122	0.09166	0.07723	0.198179	3.20033	0.058767008	0.058767008
11.1111	1087	9352	1081	0.080953	0.20163	0.07717	0.122588	0.122588	0.487742	0.02919843	0.02919843
Control	Cell 11714			Blank	0	0	0	0	0	IC50	641.7736748

33

Second batch: (r9IFN-co)

[illegible]

15

Third batch: (451FN-co)

Concentration (* 10 ⁶ /ml)	Inhibition effect to HBeAg					Accumulated inhibition rate	1- Accumulat ion	Accumulated inhibition rate
	First well	Second well	Third well	Inhibition rate First well	Second well	Average inhibition rate		
900	9702	5614	8110	0.428016	0.433204	0.52187	0.43031	1.314983
300	8914	10032	8870	0.4744723	0.40856	0.47705	0.453366	0.443859127
100	16312	12688	13934	0.398321	0.251975	0.37851	0.156271	0.402596
33.33333	15080	12814	13288	0.110954	0.244547	0.21660	0.190701	0.2463153
11.11111	21928	15366	15728	0	0.094093	0.07275	0.0055615	0.0951615
Control	Cell	17544		Blank	0	Dilution	3	IC50

Concentration (* 10 ⁶ /ml)	Inhibition effect to HBeAg					Accumulated inhibition rate	1- Accumulat ion	Accumulated inhibition rate
	First well	Second well	Third well	Inhibition rate First well	Second well	Average inhibition rate		
900	5616	6228	5346	0.496084	0.442035	0.52105	0.48651	0.761125
300	8542	8596	7096	0.234725	0.230425	0.38427	0.276474	0.2764738
100	11420	11360	11394	0	0	0	0	0
33.33333	12666	11582	13110	0	0	0	0	0
11.11111	13142	12336	13344	0	0	0	0	0
Control	Cell	11528		Blank	0	Dilution	3	IC50

HBeAg: Average IC50: 450.2434

SD: 132.315479

HBeAg: Average IC50: 649.1694

SD: 42.29580

35

Table 2: Results of inhibition rate of Icteron A (IEN-62b) to HBsAg and HBeAg

Concentration (μ M/10 ⁶ cell)	Inhibition effect to HBsAg					1- Accumulation rate	Accumulated inhibition rate
	First well	Second well	Third well	Inhibition rate First well	Inhibition rate Second well		
300	14918	11724	9950	0	0.029711	0.176529	0.068747
100	14868	16890	15182	0	0	0	0
33.3333	15780	21716	16408	0	0	0	0
11.1111	20856	15042	16108	0	0	0	0
3.703704	12083	12083	12083	0	0	0	0
Control	Cell 17544	Blank	Blank	0	0	0	0
Inhibition effect to HBeAg							
Concentration (μ M/10 ⁶ cell)	First well	Second well	Third well	Inhibition rate First well	Inhibition rate Second well	1- Accumulation rate	Accumulated inhibition rate
	First well	Second well	Third well	First well	Second well		
300	9226	8196	9658	0.152169	0.247106	0.521054	0.189295
100	10946	10340	10828	0	0.050156	0.364272	0.0184947
33.3333	12250	12390	13934	0	0	0	0
11.1111	12634	12342	12900	0	0	0	0
3.703704	10886	10886	10886	0	0	0	0
Control	Cell 10886	Blank	Blank	0	0	0	0
Inhibition effect to HBeAg							
Concentration (μ M/10 ⁶ cell)	First well	Second well	Third well	Inhibition rate First well	Inhibition rate Second well	1- Accumulation rate	Accumulated inhibition rate
	First well	Second well	Third well	First well	Second well		
300	9226	8196	9658	0.152169	0.247106	0.521054	0.189295
100	10946	10340	10828	0	0.050156	0.364272	0.0184947
33.3333	12250	12390	13934	0	0	0	0
11.1111	12634	12342	12900	0	0	0	0
3.703704	10886	10886	10886	0	0	0	0
Control	Cell 10886	Blank	Blank	0	0	0	0

Table 3: Results of inhibition rate of Infigen® to HBsAg and HBeAg
First batch: (Infigen®)

Concentration n(=10 ⁴ /ul)	First well	Second well	Third well	Inhibition effect to HBsAg			Average inhibition rate	1- Accumulation n	Accumulate d inhibition rate
				First well	Second well	Third well			
900	1417	12156	1350	0.091655	0.22068	0	0.104175	0.306157	0.1954825
300	1339	12288	1625	0.141776	0.21420	0	0.118062	0.2019822	1.777164
100	1436	18834	1419	0.079349	0	0.09024	0.056231	0.083921	2.721232
33.33333	1572	16034	1834	0	0	0	0	0.0273897	3.721232
11.11111	1750	17692	3132	0	0	0.08216	0.02739	0.02739	4.693843
Control	Cell	13602	Blank	0	0	0	0	IC50	FALSE

Concentration n(=10 ⁴ /ul)	First well	Second well	Third well	Inhibition effect to HBeAg			Average inhibition rate	1- Accumulation n	Accumulate d inhibition rate
				First well	Second well	Third well			
900	1208	11592	1223	0	0.01275	0	0.00425	0.025163	0.99575
300	1084	11484	1235	0	0.03031	0	0.010104	0.0205125	1.985646
100	1288	14696	1508	0	0	0	0	0.010808	2.985646
33.33333	1503	12928	1302	0	0	0	0	0.0108081	3.985646
11.11111	1719	11964	1150	0.004137	0	0.02828	0.01808	0.01808	4.974837
Control	Cell	11843	Blank	0	0	0	0	IC50	FALSE

Table 3: Results of inhibition rate of Interferon* to HBsAg and HBeAg
First batch: [Interferon*]

Concentration n(×10 ⁴ IU/ml)	Inhibition rate			Inhibition effect to HBsAg			Inhibition effect to HBeAg		
	First well	Second well	Third well	First well	Second well	Third well	Average inhibition rate	Accumulation n	1- 2- 3- 4- 5- 6- 7- 8- 9- 10- 11- 12- 13- 14- 15- 16- 17- 18- 19- 20- 21- 22- 23- 24- 25- 26- 27- 28- 29- 30- 31- 32- 33- 34- 35- 36- 37- 38- 39- 40- 41- 42- 43- 44- 45- 46- 47- 48- 49- 50- 51- 52- 53- 54- 55- 56- 57- 58- 59- 60- 61- 62- 63- 64- 65- 66- 67- 68- 69- 70- 71- 72- 73- 74- 75- 76- 77- 78- 79- 80- 81- 82- 83- 84- 85- 86- 87- 88- 89- 90- 91- 92- 93- 94- 95- 96- 97- 98- 99- 100- 101- 102- 103- 104- 105- 106- 107- 108- 109- 110- 111- 112- 113- 114- 115- 116- 117- 118- 119- 120- 121- 122- 123- 124- 125- 126- 127- 128- 129- 130- 131- 132- 133- 134- 135- 136- 137- 138- 139- 140- 141- 142- 143- 144- 145- 146- 147- 148- 149- 150- 151- 152- 153- 154- 155- 156- 157- 158- 159- 160- 161- 162- 163- 164- 165- 166- 167- 168- 169- 170- 171- 172- 173- 174- 175- 176- 177- 178- 179- 180- 181- 182- 183- 184- 185- 186- 187- 188- 189- 190- 191- 192- 193- 194- 195- 196- 197- 198- 199- 200- 201- 202- 203- 204- 205- 206- 207- 208- 209- 210- 211- 212- 213- 214- 215- 216- 217- 218- 219- 220- 221- 222- 223- 224- 225- 226- 227- 228- 229- 230- 231- 232- 233- 234- 235- 236- 237- 238- 239- 240- 241- 242- 243- 244- 245- 246- 247- 248- 249- 250- 251- 252- 253- 254- 255- 256- 257- 258- 259- 260- 261- 262- 263- 264- 265- 266- 267- 268- 269- 270- 271- 272- 273- 274- 275- 276- 277- 278- 279- 280- 281- 282- 283- 284- 285- 286- 287- 288- 289- 290- 291- 292- 293- 294- 295- 296- 297- 298- 299- 300- 301- 302- 303- 304- 305- 306- 307- 308- 309- 310- 311- 312- 313- 314- 315- 316- 317- 318- 319- 320- 321- 322- 323- 324- 325- 326- 327- 328- 329- 330- 331- 332- 333- 334- 335- 336- 337- 338- 339- 340- 341- 342- 343- 344- 345- 346- 347- 348- 349- 350- 351- 352- 353- 354- 355- 356- 357- 358- 359- 360- 361- 362- 363- 364- 365- 366- 367- 368- 369- 370- 371- 372- 373- 374- 375- 376- 377- 378- 379- 380- 381- 382- 383- 384- 385- 386- 387- 388- 389- 390- 391- 392- 393- 394- 395- 396- 397- 398- 399- 400- 401- 402- 403- 404- 405- 406- 407- 408- 409- 410- 411- 412- 413- 414- 415- 416- 417- 418- 419- 420- 421- 422- 423- 424- 425- 426- 427- 428- 429- 430- 431- 432- 433- 434- 435- 436- 437- 438- 439- 440- 441- 442- 443- 444- 445- 446- 447- 448- 449- 450- 451- 452- 453- 454- 455- 456- 457- 458- 459- 460- 461- 462- 463- 464- 465- 466- 467- 468- 469- 470- 471- 472- 473- 474- 475- 476- 477- 478- 479- 480- 481- 482- 483- 484- 485- 486- 487- 488- 489- 490- 491- 492- 493- 494- 495- 496- 497- 498- 499- 500- 501- 502- 503- 504- 505- 506- 507- 508- 509- 510- 511- 512- 513- 514- 515- 516- 517- 518- 519- 520- 521- 522- 523- 524- 525- 526- 527- 528- 529- 530- 531- 532- 533- 534- 535- 536- 537- 538- 539- 540- 541- 542- 543- 544- 545- 546- 547- 548- 549- 550- 551- 552- 553- 554- 555- 556- 557- 558- 559- 560- 561- 562- 563- 564- 565- 566- 567- 568- 569- 570- 571- 572- 573- 574- 575- 576- 577- 578- 579- 580- 581- 582- 583- 584- 585- 586- 587- 588- 589- 590- 591- 592- 593- 594- 595- 596- 597- 598- 599- 600- 601- 602- 603- 604- 605- 606- 607- 608- 609- 610- 611- 612- 613- 614- 615- 616- 617- 618- 619- 620- 621- 622- 623- 624- 625- 626- 627- 628- 629- 630- 631- 632- 633- 634- 635- 636- 637- 638- 639- 640- 641- 642- 643- 644- 645- 646- 647- 648- 649- 650- 651- 652- 653- 654- 655- 656- 657- 658- 659- 660- 661- 662- 663- 664- 665- 666- 667- 668- 669- 670- 671- 672- 673- 674- 675- 676- 677- 678- 679- 680- 681- 682- 683- 684- 685- 686- 687- 688- 689- 690- 691- 692- 693- 694- 695- 696- 697- 698- 699- 700- 701- 702- 703- 704- 705- 706- 707- 708- 709- 710- 711- 712- 713- 714- 715- 716- 717- 718- 719- 720- 721- 722- 723- 724- 725- 726- 727- 728- 729- 730- 731- 732- 733- 734- 735- 736- 737- 738- 739- 740- 741- 742- 743- 744- 745- 746- 747- 748- 749- 750- 751- 752- 753- 754- 755- 756- 757- 758- 759- 760- 761- 762- 763- 764- 765- 766- 767- 768- 769- 770- 771- 772- 773- 774- 775- 776- 777- 778- 779- 780- 781- 782- 783- 784- 785- 786- 787- 788- 789- 790- 791- 792- 793- 794- 795- 796- 797- 798- 799- 800- 801- 802- 803- 804- 805- 806- 807- 808- 809- 810- 811- 812- 813- 814- 815- 816- 817- 818- 819- 820- 821- 822- 823- 824- 825- 826- 827- 828- 829- 830- 831- 832- 833- 834- 835- 836- 837- 838- 839- 840- 841- 842- 843- 844- 845- 846- 847- 848- 849- 850- 851- 852- 853- 854- 855- 856- 857- 858- 859- 860- 861- 862- 863- 864- 865- 866- 867- 868- 869- 870- 871- 872- 873- 874- 875- 876- 877- 878- 879- 880- 881- 882- 883- 884- 885- 886- 887- 888- 889- 890- 891- 892- 893- 894- 895- 896- 897- 898- 899- 900- 901- 902- 903- 904- 905- 906- 907- 908- 909- 910- 911- 912- 913- 914- 915- 916- 917- 918- 919- 920- 921- 922- 923- 924- 925- 926- 927- 928- 929- 930- 931- 932- 933- 934- 935- 936- 937- 938- 939- 940- 941- 942- 943- 944- 945- 946- 947- 948- 949- 950- 951- 952- 953- 954- 955- 956- 957- 958- 959- 960- 961- 962- 963- 964- 965- 966- 967- 968- 969- 970- 971- 972- 973- 974- 975- 976- 977- 978- 979- 980- 981- 982- 983- 984- 985- 986- 987- 988- 989- 990- 991- 992- 993- 994- 995- 996- 997- 998- 999- 1000- 1001- 1002- 1003- 1004- 1005- 1006- 1007- 1008- 1009- 1010- 1011- 1012- 1013- 1014- 1015- 1016- 1017- 1018- 1019- 1020- 1021- 1022- 1023- 1024- 1025- 1026- 1027- 1028- 1029- 1030- 1031- 1032- 1033- 1034- 1035- 1036- 1037- 1038- 1039- 1040- 1041- 1042- 1043- 1044- 1045- 1046- 1047- 1048- 1049- 1050- 1051- 1052- 1053- 1054- 1055- 1056- 1057- 1058- 1059- 1060- 1061- 1062- 1063- 1064- 1065- 1066- 1067- 1068- 1069- 1070- 1071- 1072- 1073- 1074- 1075- 1076- 1077- 1078- 1079- 1080- 1081- 1082- 1083- 1084- 1085- 1086- 1087- 1088- 1089- 1090- 1091- 1092- 1093- 1094- 1095- 1096- 1097- 1098- 1099- 1100- 1101- 1102- 1103- 1104- 1105- 1106- 1107- 1108- 1109- 1110- 1111- 1112- 1113- 1114- 1115- 1116- 1117- 1118- 1119- 1120- 1121- 1122- 1123- 1124- 1125- 1126- 1127- 1128- 1129- 1130- 1131- 1132- 1133- 1134- 1135- 1136- 1137- 1138- 1139- 1140- 1141- 1142- 1143- 1144- 1145- 1146- 1147- 1148- 1149- 1150- 1151- 1152- 1153- 1154- 1155- 1156- 1157- 1158- 1159- 1160- 1161- 1162- 1163- 1164- 1165- 1166- 1167- 1168- 1169- 1170- 1171- 1172- 1173- 1174- 1175- 1176- 1177- 1178- 1179- 1180- 1181- 1182- 1183- 1184- 1185- 1186- 1187- 1188- 1189- 1190- 1191- 1192- 1193- 1194- 1195- 1196- 1197- 1198- 1199- 1200- 1201- 1202- 1203- 1204- 1205- 1206- 1207- 1208- 1209- 1210- 1211- 1212- 1213- 1214- 1215- 1216- 1217- 1218- 1219- 1220- 1221- 1222- 1223- 1224- 1225- 1226- 1227- 1228- 1229- 1230- 1231- 1232- 1233- 1234- 1235- 1236- 1237- 1238- 1239- 1240- 1241- 1242- 1243- 1244- 1245- 1246- 1247- 1248- 1249- 1250- 1251- 1252- 1253- 1254- 1255- 1256- 1257- 1258- 1259- 1260- 1261- 1262- 1263- 1264- 1265- 1266- 1267- 1268- 1269- 1270- 1271- 1272- 1273- 1274- 1275- 1276- 1277- 1278- 1279- 1280- 1281- 1282- 1283- 1284- 1285- 1286- 1287- 1288- 1289- 1290- 1291- 1292- 1293- 1294- 1295- 1296- 1297- 1298- 1299- 1300- 1301- 1302- 1303- 1304- 1305- 1306- 1307- 1308- 1309- 1310- 1311- 1312- 1313- 1314- 1315- 1316- 1317- 1318- 1319- 1320- 1321- 1322- 1323- 1324- 1325- 1326- 1327- 1328- 1329- 1330- 1331- 1332- 1333- 1334- 1335- 1336- 1337- 1338- 1339- 1340- 1341- 1342- 1343- 1344- 1345- 1346- 1347- 1348- 1349- 1350- 1351- 1352- 1353- 1354- 1355- 1356- 1357- 1358- 1359- 1360- 1361- 1362- 1363- 1364- 1365- 1366- 1367- 1368- 1369- 1370- 1371- 1372- 1373- 1374- 1375- 1376- 1377- 1378- 1379- 1380- 1381- 1382- 1383- 1384- 1385- 1386- 1387- 1388- 1389- 1390- 1391- 1392- 1393- 1394- 1395- 1396- 1397- 1398- 1399- 1400- 1401- 1402- 1403- 1404- 1405- 1406- 1407- 1408- 1409- 1410- 1411- 1412- 1413- 1414- 1415- 1416- 1417- 1418- 1419- 1420- 1421- 1422- 1423- 1424- 1425- 1426- 1427- 1428- 1429- 1430- 1431- 1432- 1433- 1434- 1435- 1436- 1437- 1438- 1439- 1440- 1441- 1442- 1443- 1444- 1445- 1446- 1447- 1448- 1449- 1450- 1451- 1452- 1453- 1454- 1455- 1456- 1457- 1458- 1459- 1460- 1461- 1462- 1463- 1464- 1465- 1466- 1467- 1468- 1469- 1470- 1471- 1472- 1473- 1474- 1475- 1476- 1477- 1478- 1479- 1480- 1481- 1482- 1483- 1484- 1485- 1486- 1487- 1488- 1489- 1490- 1491- 1492- 1493- 1494- 1495- 1496- 1497- 1498- 1499- 1500- 1501- 1502- 1503- 1504- 1505- 1506- 1507- 1508- 1509- 1510- 1511- 1512- 1513- 1514- 1515- 1516- 1517- 1518- 1519- 1520- 1521- 1522- 1523- 1524- 1525- 1526- 1527- 1528- 1529- 1530- 1531- 1532- 1533- 1534- 1535- 1536- 1537- 1538- 1539- 1540- 1541- 1542- 1543- 1544- 1545- 1546- 1547- 1548- 1549- 1550- 1551- 1552- 1553- 1554- 1555- 1556- 1557- 1558- 1559- 1560- 1561- 1562- 1563- 1564- 1565- 1566- 1567- 1568- 1569- 1570- 1571- 1572- 1573- 1574- 1575- 1576- 1577- 1578- 1579- 1580- 1581- 1582- 1583- 1584- 1585- 1586- 1587- 1588- 1589- 1590- 1591- 1592- 1593- 1594- 1595- 1596- 1597- 1598- 1599- 1600- 1601- 1602- 1603- 1604- 1605- 1606- 1607- 1608- 1609- 1610- 1611- 1612- 1613- 1614- 1615- 1616- 1617- 1618- 1619- 1620- 1621- 1622- 1623- 1624-<

Second batch: (Infergen[®])

Inhibition effect to Hsag									
Concentratio n (x10 ¹⁰ /ml)	First well	Second well	Third well	Inhibition rate			Average inhibition rate	Accumulation n	Accumulate d inhibition rate
				First well	Second well	Third well			
900	6278	6376	6408	0.200951	0.18756	0.18348	0.190367	0.274635	0.809633
300	7692	9092	6394	0.613877	0	0.18527	0.563873	0.0842678	0.0618100
100	8960	7474	8190	0	0.04765	0	0.015885	0.015885	0.00579485
33.33333	8536	8144	8682	0	0	0	0	0	0
11.11111	7048	7848	7848	0	0	0	0	0	0
Control	Cell	7848		Blank	0	0	Dilution	3	IC50
									FALSE

Inhibition effect to Hsag									
Concentratio n (x10 ¹⁰ /ml)	First well	Second well	Third well	Inhibition rate			Average inhibition rate	Accumulation n	Accumulate d inhibition rate
				First well	Second well	Third well			
900	1236	1268	1227	0.036471	0.0385	0.04318	0.04004	0.140162	0.958996
300	1159	12708	1371	0.096507	0.00935	0	0.035267	0.0991561	0.0490186
100	1244	13468	1398	0.029623	0	0	0.008674	0.063671	0.0214954
33.33333	1261	11346	1244	0.016526	0.11552	0.02293	0.053986	0.053986	0.01379630
11.11111	1282	12828	1282	0	0	0	0	0	0
Control	Cell	12828		Blank	0	0	Dilution	3	IC50
									FALSE

38

Third batch: [interferon[®]]

Concentration (x10 ⁶ /ml)	Inhibition effect to HBsAg			Inhibition effect to HBsAg			Average inhibition rate	Accumulation	3- Accumulation	Accumulation rate
	First well	Second well	Third well	First well	Second well	Third well				
900	7240	6642	6158	0.06455	0.1413	0.20433	0.136951	0.217399	0.85049	0.201217
300	11072	8786	7592	0	0	0.10526	0.03609	0.0804479	1.82696	0.0121765
300	7016	9726	7552	0.09354	0	0.02428	0.039276	0.044358	2.707693	0.015530
33.33333	7622	8866	8676	0.01524	0	0	0.005082	0.0050818	3.782601	0.003116
11.11111	7740	7740	7740	0	0	0	0	0	4.782601	0
Control	Cell	7740		Blank	0	0	Dilution	3	IC50	FALSE

Concentration (x10 ⁶ /ml)	Inhibition effect to HBsAg			Inhibition effect to HBsAg			Average inhibition rate	Accumulation	f- Accumulation	Accumulation rate
	First well	Second well	Third well	First well	Second well	Third well				
900	11048	11856	11902	0.04775	0	0	0.015917	0.015917	0.984083	0.0159167
300	13454	12896	11798	0	0	0	0	0	1.984083	0
100	12846	13160	12546	0	0	0	0	0	2.984083	0
33.33333	12680	12458	12260	0	0	0	0	0	3.984083	0
11.11111	11602	11602	11602	0	0	0	0	0	4.984083	0
Control	Cell	11602		Blank	0	0	Dilution	3	IC50	FALSE

HBsAg: Average IC50: 0 SD: 0

HBsAg: Average IC50: 0 SD: 0

39

EXAMPLE 5**Preparation of rSIFN-co****Preparation of lyophilized injection****Lyophilized powder**

Stock Solution of 34.5 µg/ml

rSIFN-co

PB (pH7.0) 10mmol/L

Glycine 0.4mol/L

5

Preparation technique: Weigh materials according to recipe. Dissolve with sterile and pyrogen-free water. Filter through 0.22µm membrane to de-bacterialize, preserve at 6-10°C. Fill in vials after affirming it is sterile and pyrogen-free, 0.3 ml /vial or 0.5 ml/vial, and lyophilize in freeze dryer.

10

Preparation of liquid injection**Solution**

Stock Solution of 34.5 µg/ml

rSIFN-co

PB (pH7.0) 26mmol/L

NaCl 0.1mol/L

15

Preparation: Weigh materials according to recipe. Add to desired level with sterile and pyrogen-free water. Filter through 0.22µm membrane to de-bacterialize, preserve at 6-10°C. Fill in airtight vial after affirming it is sterile and non-pyrogen at 0.3 ml /vial or 0.5 ml/vial. Storage at 2-10°C, and protect from light.

20

EXAMPLE 6**Acute Toxicity of rSIFN-co**

25

Treat mice with large dose (150µg/kg, equal to 1000 times of the normal dose per kilo used in treatment of adult patients) of rSIFN-co at one time by intramuscular injection. Then, observe and record their deaths and toxic

reactions. Results show that: 24 hours after injection, no abnormal reaction had been recorded. The organs of the animals which had been selected to be killed also had no signs of abnormal changes. Those remaining mice were all kept alive and were normal after two weeks. The weights of mice in the experimental group and control group all increased, and the ratio of increase had no obvious difference between the two groups ($P>0.05$) according to their weights on the fourteenth day. No abnormal changes were seen from the main organs of those mice after two weeks.

1. Experimental material

1.1 Animals

40 healthy adult mice, weighing 18-22g, half male and half female, qualified by Sichuan experiment animal control center.

2.2 Medicines

rSIFN-co (Provided by Sichuan Huiyang Life-engineering Ltd.) sterilized solution, 0.15 mg/ml, Lot: 981201
rSIFN-co was administered i.m. in saline.

2. Method

Separate the 40 mice into two groups randomly, one for experimental medicine, another for control. Inject medicines or saline at the same ratio (0.1 ml/10 g) through muscle to each mouse according to which group they belong. (150 μ g/kg of rSIFN-co for experimental group; and saline for control group). After injection, observe and record acute toxicity shown in mice. Kill half of the mice (male and female each half) to check whether there were any abnormal pathologic changes in their main organs, such as heart, spleen, liver, lung, kidney, adrenal gland, stomach, duodenum, etc. after 24 hours. Those remains were kept and observed until the fourteenth day. Weigh all mice, kill them, and then observe the appearance of the organs listed above to see if there are any abnormalities. Take

pathological tissue and examine it, using the examination to assess the difference in weight increases in the two groups.

3. Results

Results show that there was no acute toxicity seen after all mice were treated with i.m. rSIFN-co with 150 µg/kg at a time, equal to 1000 times the normal dose per kilo used in treatment of adult patients. In the 14 days after injection, all mice lived well. They ate, drank, exercised, and excreted normally and showed normal hair conditions. None of them died. The observation of the main organs of the randomly selected mice shows no abnormal changes 24 hours after injection. 14 days after injection, all remaining mice were killed. Autopsies also showed no changes. The weights of mice in the two groups all increased, but no obvious difference was shown when assessed with statistic method ($p > 0.05$). See Table 1:

Table 1 Influence to weights of mice after injection of rSIFN-co

Group	Dose	Animal	Weights before injection (g)	Weights after injection (g)	Increased value of weights (g)
Control	0	20	19.8 ± 1.7	30.8 ± 2.8	11.0 ± 2.9
rSIFN-co	150	20	19.4 ± 1.7	32.1 ± 3.3	12.7 ± 4.3

3. Conclusion

Under conditions of this experiment, there were no toxic reactions in all mice after injection of rSIFN-co with 150 µg/kg. The conclusion can be reached that the maximum tolerable dose of i.m. in mice is 150 µg/kg, which is equal to 1000 times the normal dose per kilo used in treatment of adult patients.

EXAMPLE 7

The clinic effects of recombinant super-compound interferon (rSIFN-co)

The recombinant super-compound interferon (rSIFN-co) is an invention for viral disease therapy, especially for hepatitis. Meanwhile, it can inhibit the activity of EB viruses, VSV, Herpes simplex viruses, coronaviruses, measles viruses et al. Using Wish cells /VSV system as the assay for anti-virus activity, the results showed that: the other rIFN, was 9×10^3 IU/mg, Intron A was 2.0×10^3 IU/mg and rSIFN-co was 9×10^3 IU/mg. The anti-viral activity of rSIFN-co is much higher than these of the former two.

Under the permission of the State Food and Drug Administration (SFDA), People's Republic of China, the clinical trials have taken place in Western China Hospital of Sichuan University, the Second Hospital of Chongqing Medical University and the First Affiliated Hospital of Zhejiang University School of Medicine since February 2003. The clinical treatment which focuses on the Hepatitis B is conducted under the guidance of the mutilcenter, double-blind random test. IFN- α 1b was used as control, and the primary results showed the following:

The effect of rSIFN-co compared with IFN- α 1b in the treatment of chronic active Hepatitis B

1. Standard of patients selection: The standard 1-4 are effective to both treatment with rSIFN-co (9 μ g) and IFN- α 1b (3MU, 50 μ g), and the standard 1-5 are for rSIFN-co (15 μ g) treatment.

1). Age: 18-65

2). HBsAg test positive last over six months, HBeAg test positive, PCR assay, HBV-DNA copies $\geq 10^3$ /ml

3). ALT \geq two times of the normal value

4). Never received IFN treatment; or those received the Lamivudine treatment but failed or relapsed

5) Once received other IFNs (3MU or 5MU) treatment six months ago, following the standard of SFDA but failed or relapsed

2. Evaluation of the effects:

In reference to the recommendations from the Tenth China National Committee of Virus Hepatitis and Hepatopathy, the effects were divided into three degrees according to the ALT level, HBV-DNA and HBeAg tests.

Response: ALT normal level, HBV-DNA negative, HBeAg negative

Partial response: ALT normal level, HBV-DNA or HBeAg negative

Non response: ALT, HBV-DNA and HBeAg unchanged

The response and partial response groups consider as effective cases.

3. Results of clinic trial:

Group A: treatment with rSIFN-co (9 μ g)

Group B: treatment with IFN- α 1b (5MU, 50 μ g)

Period	Group	Medicine	Cases	Effective Rate	HBeAg Transfer to negative rate	HBeAg Transfer to negative rate	HBV-DNA Transfer to negative rate	Heptal function Recover rate
8-12 weeks	A	rSIFN-co (9 μ g)	32	46.88 (15)	9.38 (3)	28.12 (9)	37.50 (12)	84.38 (27)
	B	IFN- α 1b (5MU, 50 μ g)	32	21.88 (7)	0.00 (0)	9.38 (3)	15.62 (5)	56.25 (18)
16-24 weeks	A	rSIFN-co (9 μ g)	64	54.69 (35)	7.81 (5)	25.00 (16)	34.38 (22)	90.62 (58)
	B	IFN- α 1b (5MU, 50 μ g)	64	25.00 (16)	0.00 (0)	9.38 (6)	18.75 (12)	78.13 (50)

In Group C, the cases were chronic active Hepatitis B treatment with other IFNs (3MU or 5MU) before but failed or relapsed and treated with rSIFN-co (15 μ g), subcutaneous injection, every one day, last 24 weeks. The total cases are 13. After 12 weeks treatment, 7 of 13 (53.85%) were effective. 3 of 13 (23.08%) HBeAg transferred to negative: 7 of 13 (53.85%) HBV-DNA transferred to negative: 11 of 13

(84.62%) hepatic functions recovered to normal.

4. The side effects of rSIFN-co compared with IFN- α 1b in the treatment

- 5 The side effects of IFN include fever, nausea, myalgia, anorexia, hair loss, leucopenia and thrombocytopenia, etc. The maximum dose of IFN- α 1b is 5MIU per time; the routine dose is 3 MIU. When taken the routine dose, 90% patients have I- II degree (WHO standard) side effects. They are
- 10 fever lower than 38°C, nausea, myalgia, anorexia, etc. When taken at maximum dose, the rate of side effects do not rise obviously, but are more serious. The maximum dose of rSIFN-co is 24 μ g, subcutaneous injection, every one day for 3 months. The routine dose is .9 μ g. When routine doses were
- 15 used, less than 50% patients have I-II degree (WHO standard) side effects, including fever below 38°C, nausea, myalgia, anorexia, leucopenia and thrombocytopenia slightly. With maximum dosage, about 50% patients suffered from leucopenia and thrombocytopenia after using rSIFN-co one month, but
- 20 those side effects would disappear after stopping treatment for one week. It is safe for continuous use.

The observations of rSIFN-co treat Hepatitis C

1. Standard of patient's selection

- 25 1) age: 18-65
2) HCV antibody positive
3) ALT \geq 1.5 times of the normal value, last more than 6 months

2. Evaluation of the effects:

- 30 Referring to the standard of Infergen® for treatment of Hepatitis C and according to the ALT level and HCV-RNA test, divided the effects into three degree:
Response: ALT normal level, HCV-RNA negative
Partial response: ALT normal level, HCV-RNA unchanged
35 Non response: ALT and HCV-RNA unchanged

3. Effects in clinic

The clinical trial was done at the same time with Hepatitis B treatment. 46 cases received the treatment, 9 µg each time, subcutaneous injection, every day for 24 weeks. After
5 treatment, 26 of 46 (56.52%) have obvious effects, 12 of 46 (26.08%) HCV-RNA transferred to negative, 26 of 46 (56.52%) hepatic functions recovered to normal.

10 When used in this specification and claims, the terms "comprises" and "comprising" and variations thereof mean that the specified features, steps or integers are included. The terms are not to be interpreted to exclude the presence of other features, steps or components.

What is claimed is:

1. A recombinant super-compound interferon or a functional equivalent thereof with changed spatial configuration and improved efficacy.
2. The interferon of claim 1, wherein the interferon is either α , β , or ω .
3. The interferon of claim 1, wherein the interferon has higher efficacy than the interferon described in U.S. Patent Nos. 4,695,623 or 4,897,471, or with unique secondary or tertiary structure.
4. The super-compound interferon of claim 1, wherein the spatial configuration change is the result of changes of its production process.
5. A super-compound interferon of claim 1, produced by a high efficiency expression system which uses a special promoter, or the promoter is P_{lac}.
6. The super-compound interferon of claim 4, wherein its gene is artificially synthesized cDNA with adjustment of its sequence from the wild-type according to codon preference of *E. coli*.
7. The super-compound interferon of claim 1, which possesses anti-viral or anti-tumor activity.
8. The super-compound interferon of claim 7, wherein the virus diseases comprises Hepatitis A, Hepatitis B, Hepatitis C, other types of hepatitis, infections caused by Epstein-Barr virus, Cytomegalovirus, herpes simplex viruses, other herpes viruses, papovaviruses, poxviruses, picornaviruses, adenoviruses, rhinoviruses, human T cell leukemia viruses I, human

T cell leukaemia viruses II, or human T cell leukemia viruses III.

9. The super-compound interferon of claim 8, which
5 directly inhibits the DNA duplication and secretion of
HBsAg and HBeAg of Hepatitis B Virus.

10. An artificial gene codes for the super-compound
10 interferon or its equivalent of claim 1.

11. A vector comprising the gene of claim 10, an
expression system comprising the said vector, or a
host cell comprising the said vector.

12. A process for production of recombinant super-
compound interferon comprising introducing an
artificial gene with selected codon preference into an
appropriate host, culturing said introduced host in an
20 appropriate condition for the expression of said
compound interferon and harvesting the expressed
compound interferon.

13. The process for production of claim 12, comprising
25 extraction of super-compound interferon from
fermentation broth, collection of inclusion body,
denaturation and renaturation of the harvested
protein, or wherein the process maintains the high
efficacy even when the super-compound interferon is
used with an agent and in a particular concentration,
or comprising separation and purification of the
super-compound interferon, or comprising
30 lyophilization of the purified super-compound
interferon, or comprising production of liquid
injection of super-compound interferon.

14. The produced super-compound interferon by the process

of any of the claims 12-13.

15. A composition comprising the recombinant super-compound interferon of claim 1 and a suitable carrier.

16. A pharmaceutical composition comprising the recombinant super-compound interferon of claim 1 and a pharmaceutically acceptable carrier.

17. A method for treating viral diseases or tumor in a subject comprising administering to the subject an effective amount of the super-compound interferon of claim 1.

18. The method of claim 17 wherein the viral diseases is Hepatitis A, Hepatitis B, Hepatitis C, other types of hepatitis, infections of viruses caused by Epstein-Barr virus, Cytomegalovirus, herpes simplex viruses, or other type of herpes viruses, papovaviruses, poxviruses, picornaviruses, adenoviruses, rhinoviruses, human T cell leukaemia viruses I, or human T cell leukaemia viruses II, or human T cell leukemia virus III.

19. The method of claim 17 wherein super-compound interferon was administered via oral, vein injection, muscle injection, peritoneal injection, subcutaneous injection, nasal, mucosal administration, by inhalation via an inspirator.

20. The method of claim 17 wherein super-compound interferon was administered following the protocol of injection 9 µg or 15 µg per day, 3 times a week, total 24 weeks.

Figure 1

5' 11 21 31 41 51
 +1 M G D L P G T H S L S N R R A L I L L A
 1 ATGTGCGACG TCGGCGCGAC CGATCGCTG GATAAGCGTG GTGCTGTGAT CTTGCTGAGCT
 TAGAGCGCTG AGCGGCTGCTG GGTGAGGCGAC CGATTGCGAG CAGCGABACTA GAGAGACGGA

5' 71 81 91 101 111
 +1 Q M R R I S P F S G L K D R H D F G F P
 81 CAGATGCGTC GTATGTGCGG GTTGTGCTG GCGAAGACG GTGAGGACTT CGTTTTCGCG
 GTGTAGCGAG CATAAGGCGG CAGAGAGCG GACTTTGTG CGTGTGTA GCGAAGGCG

5' 131 141 151 161 171
 +1 Q E E F D G N Q F Q K A D A I S V L N E
 121 CAGGAGAAAT TCGAGCGTA CGGTTGCGAG AAGGCTGAG GTATGCTGCT TGTGACGAA
 GTGTTGTTTA AGCTGCGATT GGTGAGGCT TGTGAGCTG GATAGAGCG AGAGTGTGCT

5' 191 201 211 221 231
 +1 M I Q G T F Y L F S T K D S S A A W D E
 181 ATGATCGACG AGAGCTGCAI CTTGTGTGCG ACCAGAGAGT CGTGTGCTG TGTGACGAA
 TACTAGTGTG TGTGAGATT GAGCAGAGCG TGCTTGTGCA GAGAGCGAG AGCGCTGTGCT

5' 251 261 271 281 291
 +1 S L L E K F Y T E L Y G G L N D L E A C
 241 TCGTGTGCG AAAATTCTA CAGGAGCTG TAGCGAGCG TCGAGAGCT GAGAGCTTGC
 AGGAGCGACG TTTTAAAGAT GTGCTGTGAG ATGTGTGCG AGTGTGTGCA GTTGTGAGCG

5' 311 321 331 341 351
 +1 V I Q E V G Y E E T P L N N V D S I L A
 301 GTTATCGAGG AGTTGTGTGT TGAGAGAACG CGGTGATGA AGTTGTGACT GATGTGTGCT
 GATAGAGTGTG TTAGCGACA ACTTGTGTG GAGGACTACT TCGAGTGTG GTAGAGCGAG

5' 371 381 391 401 411
 +1 Y K K Y F Q R I T L Y L T E K K Y S P C
 361 GTTAAAGAAAT AGTGTGAGCG TATGAGCGTG TAGCTGAGCG AAAAAAATA GTGCGGCTGCG
 CAAATTTTTA TGAGGTGTGCG ATAGTGGAGC ATGAGCTGCG TTTTITTTAT GAGGAGCGAG

5' 431 441 451 461 471
 +1 A N E V V R A E I N R S F S L S T N L O
 421 GCTTGGGAGG TGTGTGTGCG TGAAATGAGC CGTTGCTTCT CGCTGTGCG CAGAGTGTGAG
 CGAGCGCTTC AGAGAGCGAI AGTTTATGAG GAGAGGAGA GAGAGAGTGT GTTGTGAGCTG

5' 491 501
 +1 E R L R R K E
 481 GAGGTGTGCG GTGTAGAGA ATAA
 GTTGTGAGCG CAGCATTTGT TATT

Figure 2

```

5'          11      21      31      41      51
+1 M C D L P Q T H S L G N R R A L I L L A
1 ATGCGGATT TACCTGAAAC TATTGCGTT GCGACCGTC GCGCTCGGT TCTGCTGGA
TACGACTAA ATGGGTTTG AGAAGGGA CATTGGGCG GCGGAGCTA AGGAGCGGT

5'          71      81      91      1      11
+1 Q M R R I S P F S C L K D R H D F G F P
81 CAGTGGGCT GCAATTGCG GATTACGTC CCGAAGGCG GCGGAGCT CCGCTTGG
GCTAGCGAG CATAAGGCG CAAATGCG GCTTTGCG CAGCTGTA GCGGAAAGG

5'          31      41      51      61      71
+1 Q E E F D G N Q F Q K A Q A I S V L H E
121 CAGGAGGCT TCGATGGGA CATTGCG AGACCTGCG CATTCTGT ACTGACGGA
GTTCTGTA AGCTAGGCT GGTAGGTC TTGTAGTCT GTTAGGGA TGAGTGGT

5'          91      1      11      21      31
+1 M I Q Q T F N L F S T K D S S A A W D E
181 ATGATCGAC AGCGCTTA CCGTTTTC ACTAAGGA GCTGCTGC TTGCGGGA
TACTAGGTC TGTGAGCT GCGAAGGCT TATTCTGT CAGAGCGAG AACCTGCT

5'          51      61      71      81      91
+1 S L L E K F Y T E L Y Q Q L N D L E A C
241 ACCTCTGCG AGAGCTTA CAGGAGCT TATCGAGC TGAGGCT CAGAGATCC
TGGAGGCG TCTTAGGT GCGCTTGC ATAGCTGCG ACTGCTGA CATTGAGC

5'          11      21      31      41      51
+1 V I Q E V G V E E T P L M N V D S I L A
301 GCAATGCG AGCTGCTGT AGAGGCT CCGTGTA AGCTGCT TATTGGA
CATTAGGCG TTAGGGA TCTCTGTA GCGCTAC TGAGCTGAG ATAGGCTG

```

Figure 3

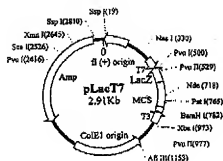


Figure 4

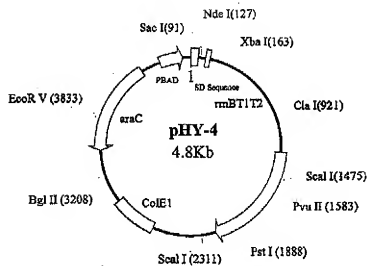


Figure 5

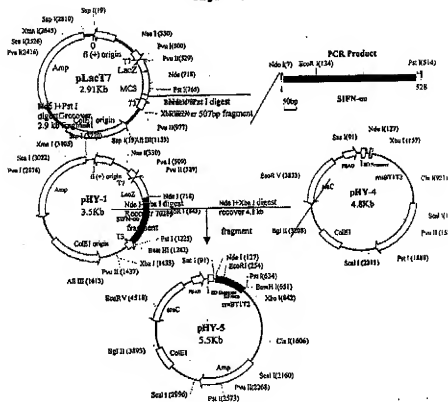


Figure 6-A

Circular Dichroism spectra

Tested by Analysis and Measurement Center of Sichuan University.

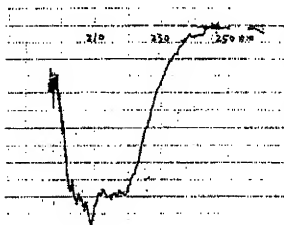


Fig 6-A Circular
Dichroism spectrum
of Inergen
Spectrum range:
250nm - 190nm
Sensitivity: 2 m°/cm
Light path: 0.20 cm
Equipment: Circular
Dichroism J-5000
Samples :contain
30µg/ml IPN-ctrl,
5.9 mg/ml of NaCl
and 3.8 mg/ml of
Na₂PO₄, pH7.0.

Figure 6-B

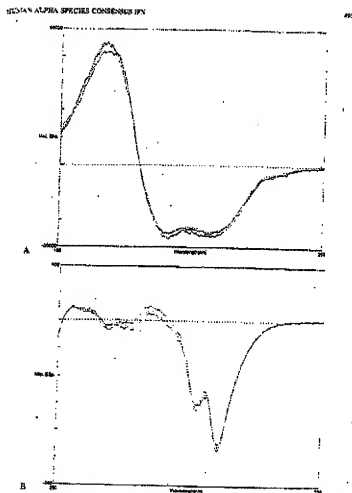


FIG. 3. Circular dichroism spectra of ribonuclease A. Ribonuclease A was purified using an ion exchange column. It shows in Figure 3. Samples were dialyzed into 10 mM sodium phosphate, pH 7.4. Measurements were made on a Jasco J-150 spectropolarimeter. In a cell thermostat at 15°C. (—) native, (---) heated (100°C), (---) heated (100°C) and cooled (15°C). A: Far UV spectra. B: Near UV spectra.

Fig 6-B Circular Dichroism spectrum of ribonuclease A [Journal of Interferon and cytokine Research. 16:489-499(1996)]

Figure 6-C

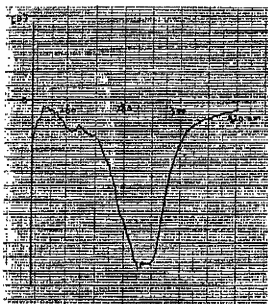


Fig 6-C Circular Dichroism spectrum of rSIFN-co

Spectrum range: 220nm-250nm

Sensitivity: 2 mdeg

Light path: 2cm

Equipment: Circular Dichroism J-500C

Samples : contain 0.5mg/ml rSIFN-co, 5.9 mg/ml of NaCl and 3.8 mg/ml of Na_2PO_4 , pH7.0.

Figure 6-D

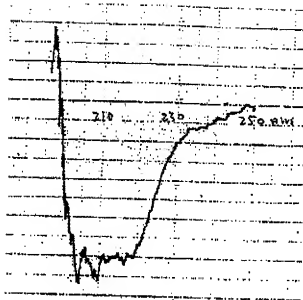


Fig 6-D Circular Dichroism spectrum of rSIFN-co

Spectrum range: 250nm - 190nm

Sensitivity: 2 m°/cm

Light path: 0.20 cm

Equipment: Circular Dichroism J-500C

Samples contain 30µg/ml rSIFN-co, 5.9 mg/ml of NaCl and 3.8 mg/ml of Na₂PO₄, pH7.0.